


Application Instructions

The G-Band identification wristband is designed for comfort and ease of use. The unique adhesive fitting and closure system makes it easy to fit the G-Band to any wrist size.

STEP 1
*Remove backing
from fitting tab (1)*

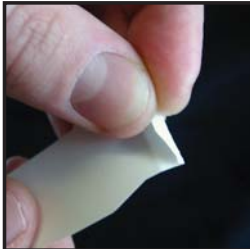

STEP 2
*Place the G-Band
around wrist*




STEP 3
*Affix the fitting
tab (1)*



STEP 4
*Remove backing
from closure tab (2)*



STEP 5
*Affix closure
tab (2)*

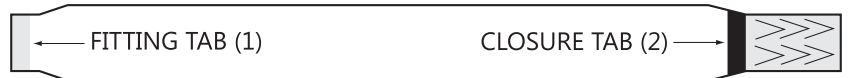

STEP 1: Remove the backing from the fitting tab (1)

Before applying to wrist, first remove the backing from the small fitting tab (1) on the back of the G-Band. Discard the backing.

STEP 2: Place the G-Band around the wrist

Place the G-Band around the wrist, with the fitting tab (1) at the bottom, facing towards you.

STEP 3: Affix the fitting tab (1)

Bring the bottom of the G-Band up around the wrist, and affix the fitting tab to the underside of the G-Band at a location where the G-Band fits correctly and comfortably around the wrist.

STEP 4: Remove the backing from the closure tab (2)

Remove the backing from the larger closure tab (2) on the back of the G-Band. Discard the backing.

STEP 5: Affix the closure tab (2)

Wrap the top of the G-Band around the wrist, and affix the closure tab (2) to the outside of the G-Band. For user comfort, make sure the closure tab is centered on the surface of the G-Band with no exposed adhesive on either side. Press on the closure tab to assure positive adhesion.

After the closure tab has been affixed, do not attempt to re-open it. The closure tab is equipped with tamper-evident slits that are designed to tear if the wristband or closure has been opened or tampered with.